




Pure-Well® Self-Sealing Mat

Patented Pattern Adhesive Sealing Mat Reseals After Penetration

Designed for square or round 96 well plates, our sealing mat is manufactured from 30 mil medical grade silicone. The tinted adhesive is patterned so that there is no adhesion exposed over the well. The adhesion also creates a firm bond that avoids cross-contamination and evaporation.

The sealing mat has superior chemical compatibility and can resist the most extreme solvents. The mat is scored for easy penetration and self-seals after needle or pipette removal.

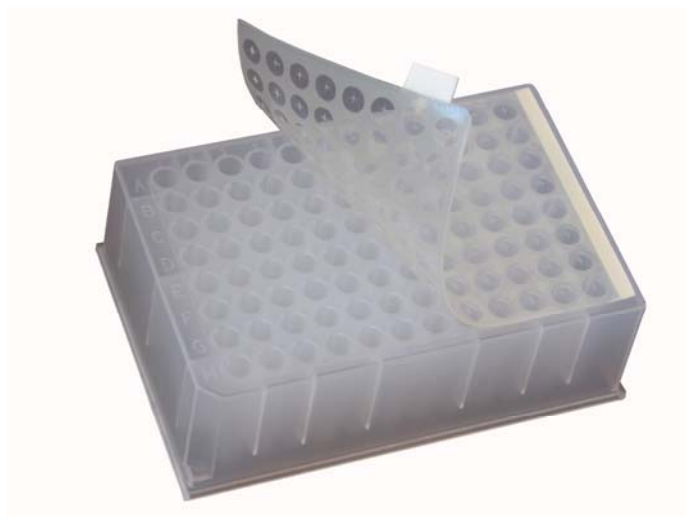

- Replaces Difficult-To-Use Dimpled Mats
- Scored for Easy Penetration
- Self Seals after Dispense or Aspiration
- Resists Clogging, Coring and Tearing
- DMSO and Solvent Resistant
- Temperature Resistant: -70°C to 140°C

Cat #	Description	QTY/PK
ES3279	Pure-Well Self Sealing Mat, 96-Well Round	10
ES3280	Pure-Well Self Sealing Mat, 96-Well Square	10