

15mL & 50mL Centrifuge Tubes



- Prevents breakage & leakage during high-speed centrifugation
- Easy-to-read graduations for accurate measurements
- Patented, extra large, smear-resistant writing panel



Centrifuge tubes offer the highest standard quality and performance.

Clarity Made from an ultra-clear, non-cytotoxic medical grade resin that allows you to see sample volume and color easily. These tubes feature highly visible printed graduations for easy filling and sample preparation.

Durability Engineered with thick, durable walls for ultra high speed centrifugation.

Reliability Patented, extra large writing panel is the largest on the market and transparent to easily view your sample through your label. This writing panel is also smear-resistant for permanent coding and won't wear off when inserting and removing your tube from the centrifuge rotor.

Security Leak resistant flat caps feature a molded-in sealing ring that keeps your closed tubes completely secure. These caps are ideal for vortexing chemicals and long rocker incubations.

Sustainability Designed for maximum product recyclability and manufactured in a solar powered facility with unique waste-free molding machines. Tube graphics and graduations are printed with UV cured links which significantly reduces the emissions of harmful volatile organic compounds (VOCs).




Special Features:




Leak resistant flat caps feature an elastomeric molded-in sealing ring that keeps your sample secure



Patented, extra large writing panel shows contents in tube clearly, labels easily, & is smear resistant for permanent coding

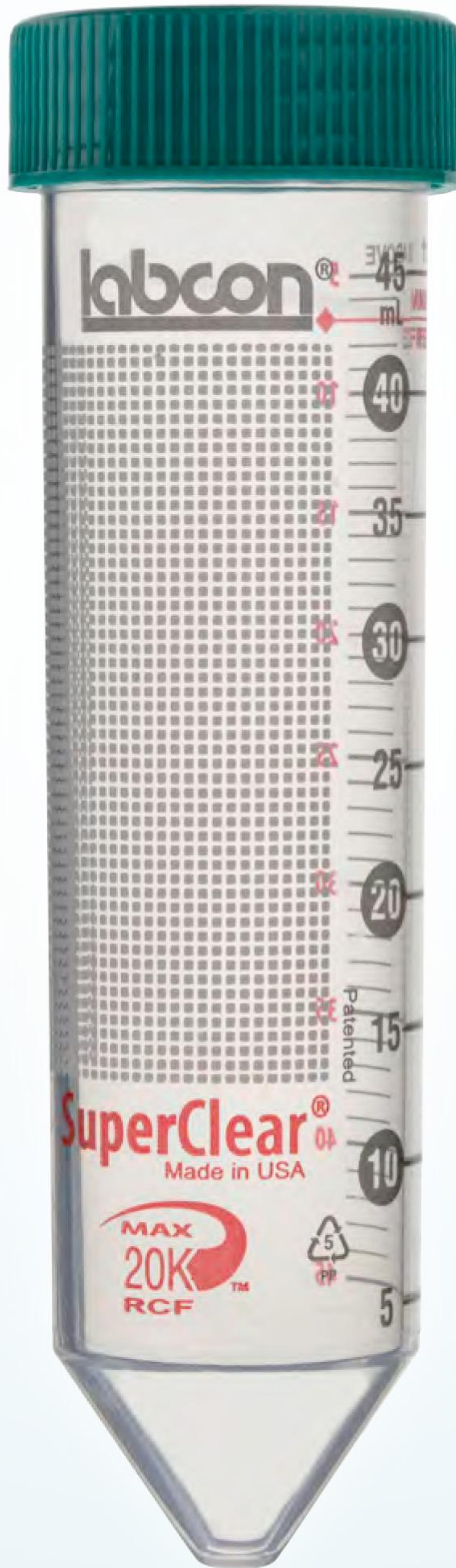

Max Freeze indicates the maximum liquid volume for samples that will be frozen

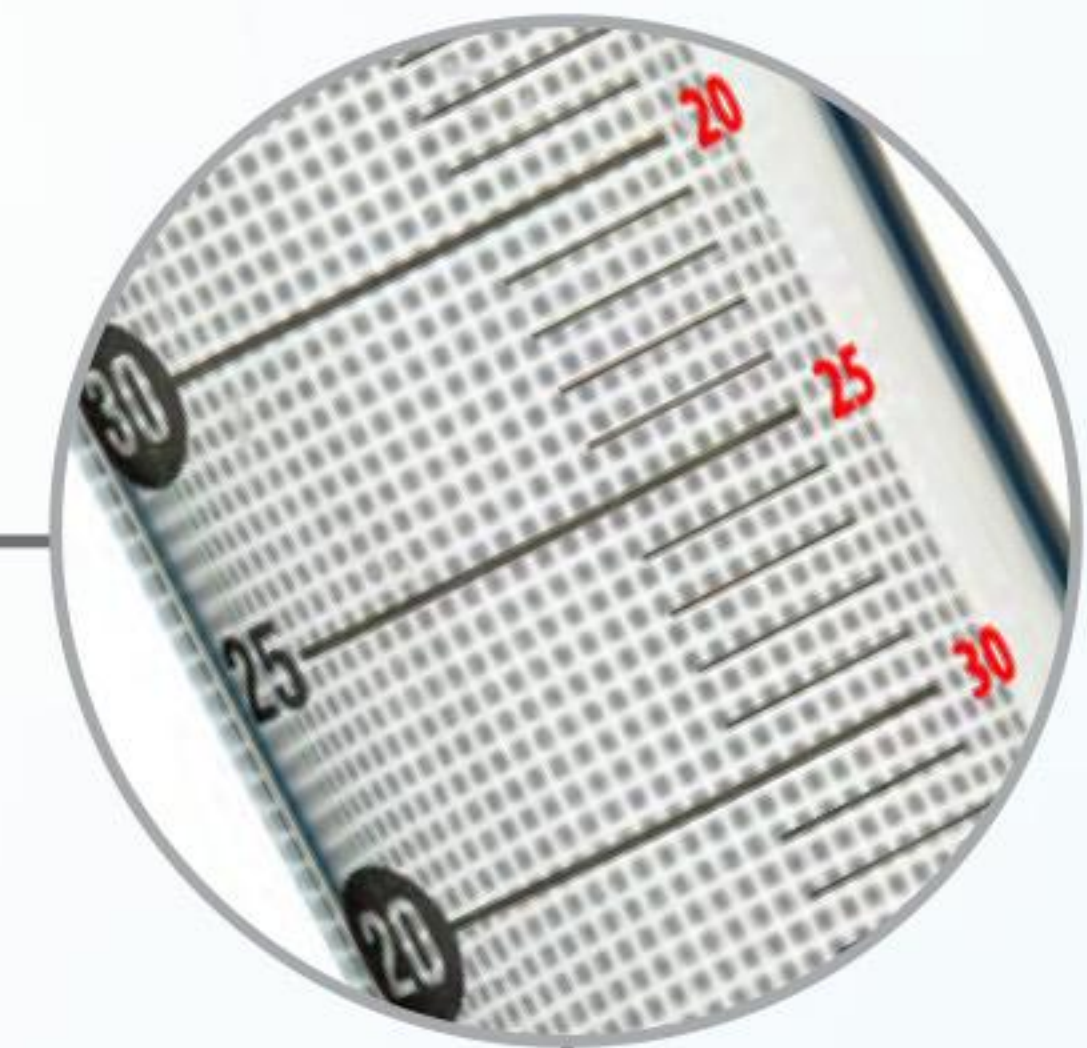

Easy-to-read graduations & reverse scale for accurate measurements



Tubes are recyclable and are offered in reusable & recyclable racks



Thick, durable, medical-grade resin withstands breakage and allows these tubes to spin at g-forces up to 233% greater than other tubes




Flat Cap Details: Cap Options

Caps feature precise vertical ridges for “easy grip” ergonomic handling, even with gloves. The low profile cap design allows for tubes to be easily stored in racks without overcrowding. Caps are autoclavable, as long as they are attached to tubes.



FLAT CAP: offers unsurpassed sealing. Both 15mL and 50mL flat caps include an integral molded-in elastomeric sealing ring that eliminates leaking. With their smooth surface, flat style caps can be easily labeled and are designed to produce an evaporation proof seal with a $\frac{3}{4}$ turn for one-handed operation plus a final $\frac{1}{4}$ turn to seal them liquid tight.

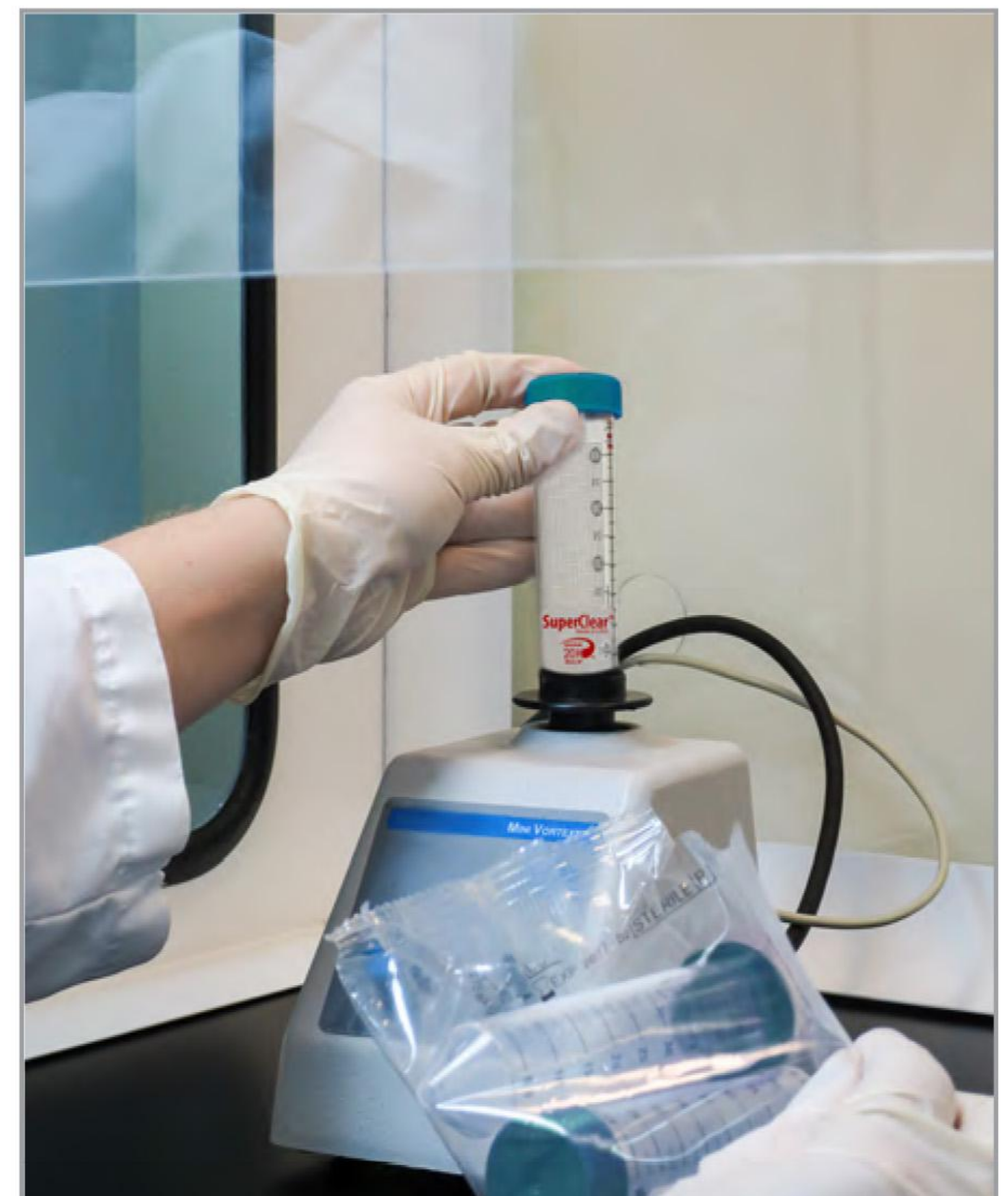

Flat caps are leak-resistant and keep your sample secure while using a vortexer or rocker.



Cat. No.	Description	UOM
ES84127	15mL Tubes, Bulk, Flat Caps, Sterile	50/pack, 10 packs/case
ES84127R	15mL Tubes, Racked, Flat Caps, Sterile	50/pack, 10 packs/case
ES84141	50mL Tubes, Bulk, Flat Caps, Sterile	50/pack, 10 packs/case
ES84141R	50mL Tubes, Racked, Flat Caps, Sterile	50/pack, 10 packs/case